


Basics for creating media

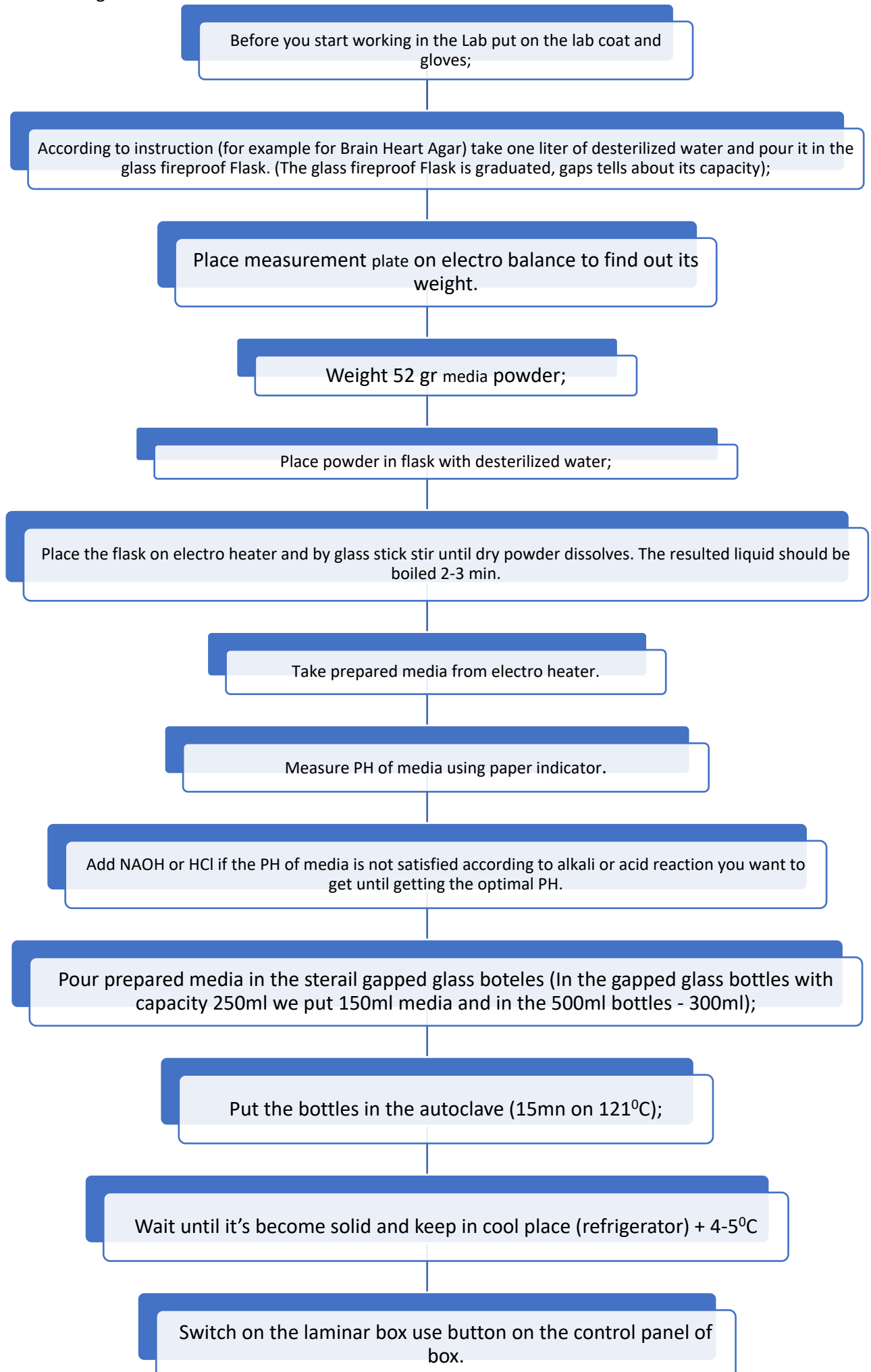

Basics for creating media

Put down the front glass.

Switch on the UV light for 30 min use another button on the control panel of box.

Afterwards switch off UV light.

Turn on air source and the light by help of another button;

Melt the media in microwave oven;

Put sterile plastic or glass Petri dishes and melted media in the laminar box.

Carefully unpack melted media bottle from foil paper and cover.

Take bottle in left hand, with right hand take up half the Petri dishes cover and pour 30ml of media, after quickly put down the cover in initial condition; .

Do the same until bottle is empty.

Media on Petri dishes become solid and should be kept with cover down to avoid accommodation of big amount of condensate on the media surface.

Label Petri dishes by marker – type of media and date.

Keep Petri dishes with media in the refrigerator on + 4-5°C.